

SHINWA EXPRESS

Application


2016
December

VOL. 97

Analysis of sertraline hydrochloride using the ULTRON ES-OVM

Sertraline hydrochloride is a medicine used to ease depressed and anxious feelings. It is also used in treatments for panic disorders and PTSD other than depressions and anxiety disorders. Shown below is a reversed-phase analysis of sertraline hydrochloride performed using the ULTRON ES-OVM.

Sertraline hydrochloride (selective serotonin reuptake inhibitor)

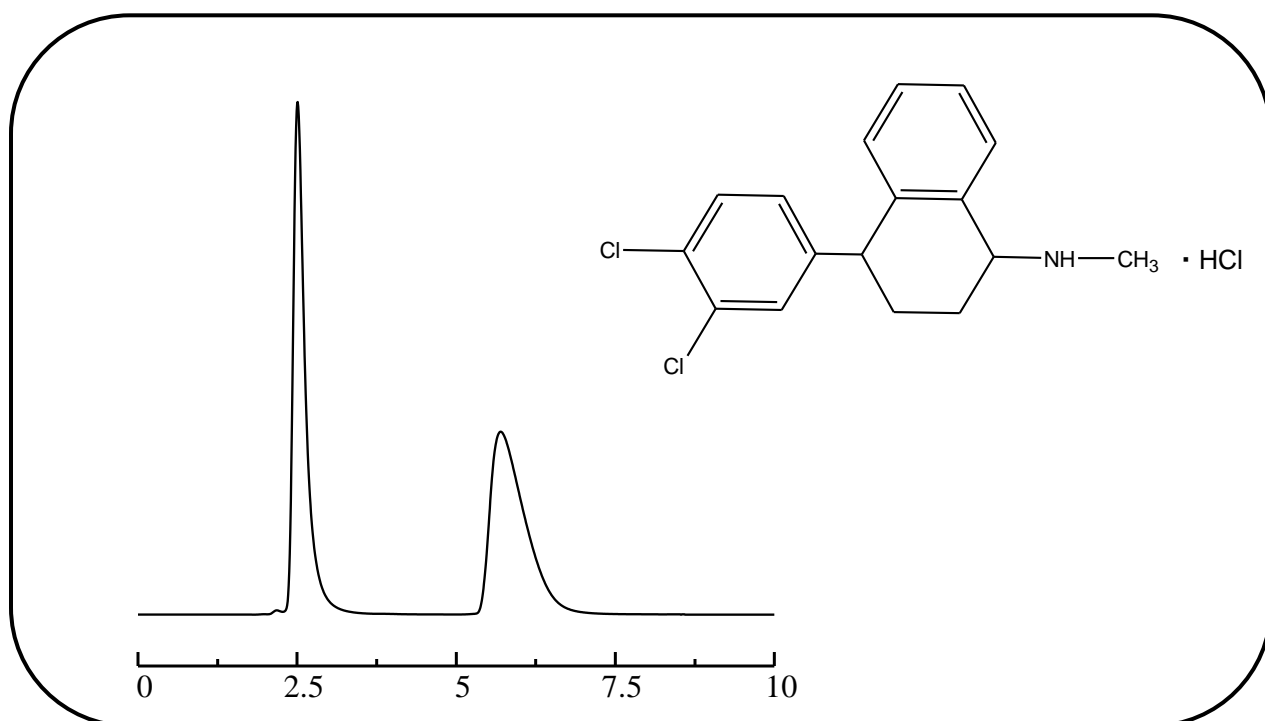

Column: ULTRON ES-OVM (5 μ m)
 Column size: 150 mm L x 4.6 mm I.D.
 Mobile phase: 20 mM KH_2PO_4 + H_3PO_4 (pH 3.0)/ CH_3CN = 85/15
 Flow rate: 1.0 mL/min

Column temp.: 25°C
 Detection: UV-220 nm
 Sample: 1000 mg/L (in methanol)
 Injection vol.: 2 μ L

●Features of the ULTRON ES-OVM

1. Particularly effective for separation of chiral compounds.
2. Employs a highly denaturation resistant protein ligand* (ovomucoid) as stationary phase.
3. Displays a wide range of chiral recognition.
4. No sample pretreatment required for chiral separations. (Blood plasma samples can be analyzed without pretreatment.)
5. Can be used for reversed-phase analyses.
6. Trace analyses possible (several nanograms).

*Ligand: A substance that binds selectively with a specific receptor.



SHINWA CHEMICAL INDUSTRIES LTD.

50-2 Kagekatsu-cho, Fushimi-ku
 Kyoto 612-8307 Japan

TEL: +81-75-621-2360 FAX: +81-75-602-2660

E-mail: info@shinwa-cpc.co.jp

Website: <http://shinwa-cpc.co.jp/en/>

■ Analytical column specifications ■

Product name	Particle size (μm)	Column dimensions Length × I.D. (mm)
ULTRON ES-OVM	5	150 × 2.0
ULTRON ES-OVM	5	150 × 4.6
ULTRON ES-OVM	5	150 × 6.0
ULTRON ES-OVM-C	5	150 × 4.6
ULTRON ES-OVM	10	250 × 4.6
ULTRON ES-OVM Prep (preparative)	10	250 × 20.0
ULTRON ES-OVM	3	Contact us

■ Guard column specifications ■

Product name	Particle size (μm)	Column dimensions Length × I.D. (mm)
ULTRON ES-OVM.G	5	10 × 4.0
ULTRON ES-OVM Prep.G (preparative guard column)	10	15 × 8.0
ULTRON ES-OVM (guard cartridges (x2))	5	5 × 2.0
Guard cartridge holder (with adaptor)	For 5 × 2.0 mm columns	
ULTRON ES-OVM (guard cartridges (x 2))	5	10 × 4.6
Guard cartridge holder (with adaptor)	For 10 × 4.6 mm columns	

● Protecting columns using guard columns

By using guard columns together with analytical columns, analytical columns can be protected from contamination, enabling longer column lifetimes. The use of guard columns is particularly recommended when analyzing samples containing large amounts of impurities.

Guard cartridges are also available when guard columns are changed frequently.

Please use guard columns with the same or smaller inner diameter than that of the analytical column being used.

● ULTRON ES-OVM-C

Clopidogrel bisulfate evaluation conditions based on USP standards employed for shipping inspections.
(Performance reports included with columns.)



Please feel free to contact us with questions related to analyses. Column screening services are also available.

Please be aware that specifications and prices are subject to change without prior notification.