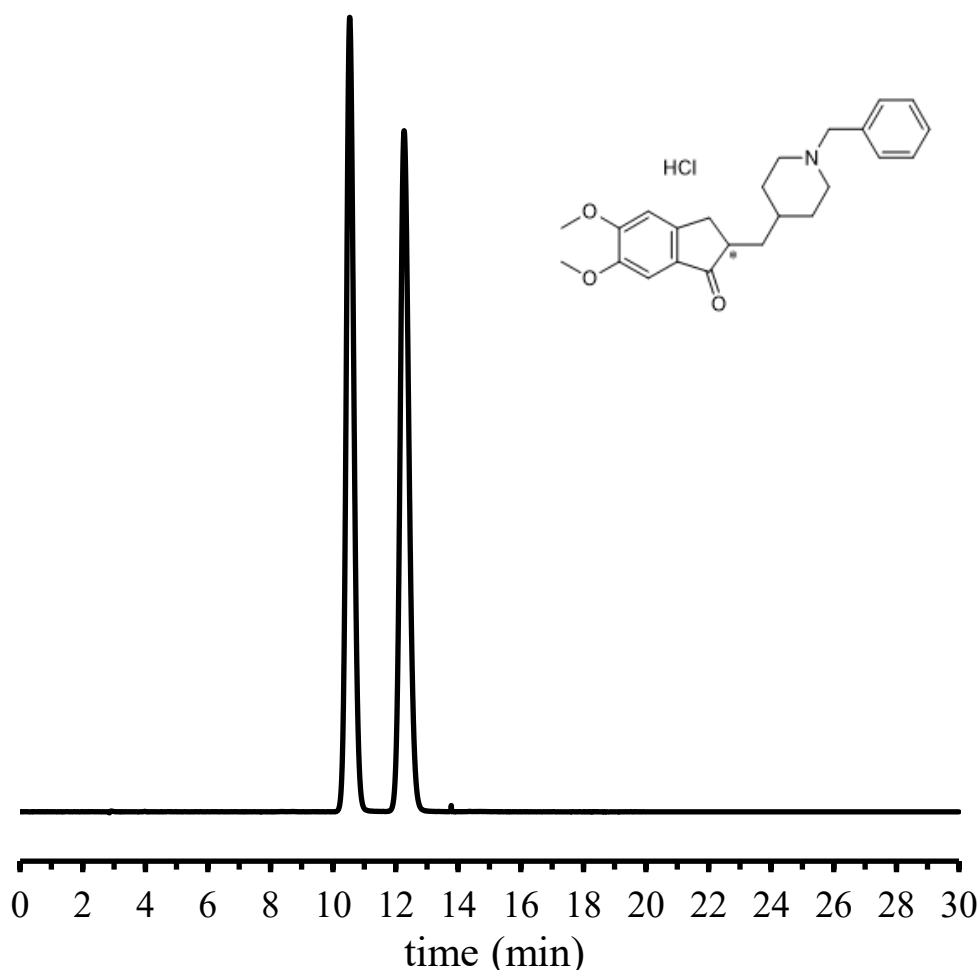


SHINWA EXPRESS

Application

VOL. 285

Analysis of Donepezil HCl using the ULTRON ES-Cel-DM

Column: ULTRON ES-Cel-DM (5 μ m)

Detection: UV-300 nm

Column size: 250 mmL x 4.6 mm I.D.

Sample: Donepezil HCl (1000 mg/L)

Mobile phase: n-Hexane/EtOH/DEA = 90/10/0.1

Diluent: mobile phase

Flow rate: 1.0 mL/min

Injection vol.: 10 μ L

Column temp.: 25°C

● Features of the ULTRON ES-Cel-DM

1. The columns are designed for the separation of various isomers
2. Corresponded to L40

**SHINWA CHEMICAL INDUSTRIES LTD.**

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>

■ Analytical column specifications ■

ULTRON ES-Amy-DM / ES-Cel-DM / ES-AmyCL-DM / ES-CelCL-DM

Product name	Particle size (μm)	Column dimensions Length x I.D. (mm)
ULTRON ES-Amy-DM	5	150×4.6
		250×4.6
	10	250×4.6
ULTRON ES-Cel-DM	5	150×4.6
		250×4.6
	10	250×4.6
ULTRON ES-AmyCL-DM	5	150×4.6
		250×4.6
	10	250×4.6
ULTRON ES-CelCL-DM	5	150×4.6
		250×4.6
	10	250×4.6

※ Please do not hesitate to contact us for the other dimensions.
Please feel free to contact us with questions related to analyses.
Column screening services are also available.



Please be aware that specifications and prices are subject to change without prior notification.