

# SHINWA EXPRESS


*Application*

2014  
J u l y  
VOL. 32

## Analysis of bupropion hydrochloride using the ULTRON ES-OVM

Bupropion hydrochloride is a dopamine reuptake inhibitor used as a medication for treating depression, and is also a nicotine antagonist used as a smoking-cessation aid in smoking cessation therapy. Shown below is a reversed-phase analysis of bupropion hydrochloride performed using the ULTRON ES-OVM.

### Bupropion hydrochloride (antidepressant, nicotine antagonist)

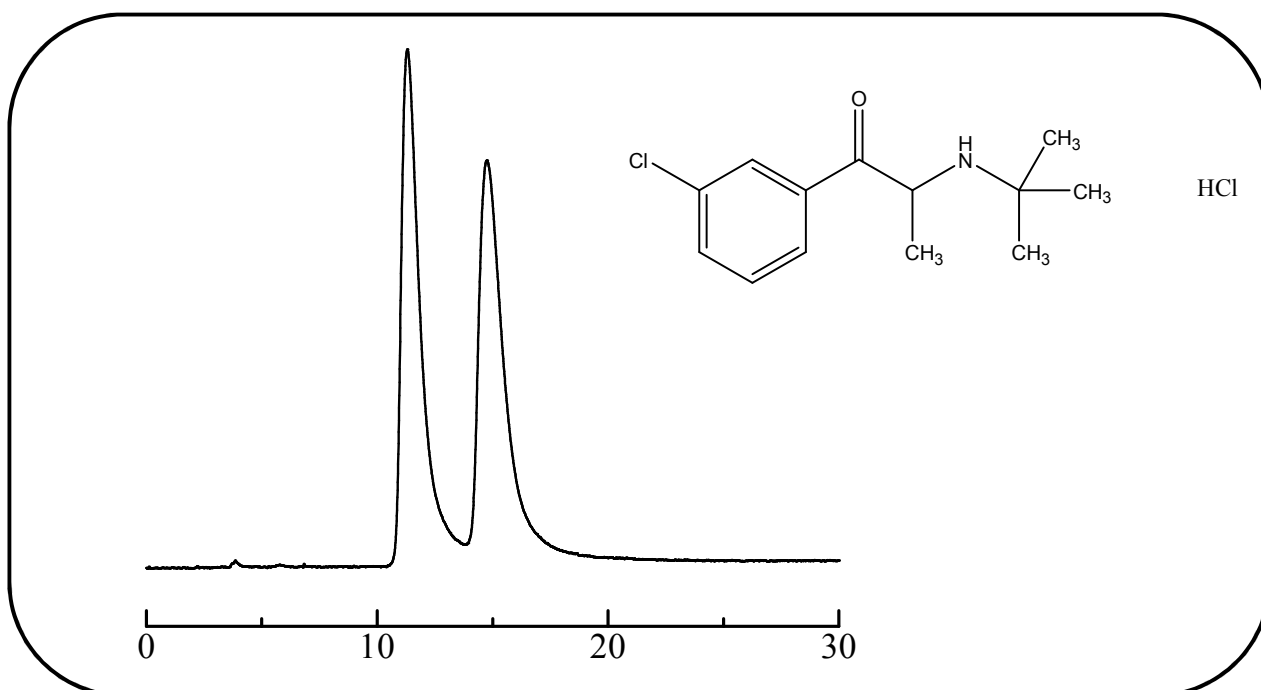

Column:	ULTRON ES-OVM (5 $\mu$ m)	Column temp.:	25°C
Column size:	150 mm L x 4.6 mm I.D.	Detection:	UV-250 nm
Mobile phase:	20 mM $\text{KH}_2\text{PO}_4$ + 20 mM $\text{K}_2\text{HPO}_4$ (pH 6.0) / $\text{C}_2\text{H}_5\text{OH}$ = 90/10	Sample:	Bupropion hydrochloride 200 mg/L in diluent (methanol : water = 1:1)
Flow rate:	1.0 mL/min	Injection vol.:	2 $\mu$ L

#### ●Features of the ULTRON ES-OVM

1. Particularly effective for separations of chiral compounds.
2. Employs a highly denaturation resistant protein ligand\* (ovomucoid) as a stationary phase.
3. Displays a wide range of chiral recognition.
4. No sample pretreatment required for chiral separations. (Blood plasma samples can be analyzed without pretreatment.)
5. Can be used for reversed-phase analyses.
6. Trace analyses possible (several nanograms).

\*Ligand: A substance that binds selectively with a specific receptor.



SHINWA CHEMICAL INDUSTRIES LTD.

50-2 Kagekatsu-cho, Fushimi-ku  
Kyoto 612-8307 Japan  
TEL: +81-75-621-2360 FAX: +81-75-602-2660  
E-mail: [info@shinwa-cpc.co.jp](mailto:info@shinwa-cpc.co.jp)  
Website: <http://shinwa-cpc.co.jp/en/>

## ■ Analytical column specifications ■

Product name	Particle size (μm)	Column dimensions Length × I.D. (mm)
ULTRON ES-OVM	5	150 × 2.0
ULTRON ES-OVM	5	150 × 4.6
ULTRON ES-OVM	5	150 × 6.0
ULTRON ES-OVM-C	5	150 × 4.6
ULTRON ES-OVM	10	250 × 4.6
ULTRON ES-OVM Prep (preparative)	10	250 × 20.0

## ■ Guard column specifications ■

Product name	Particle size (μm)	Column dimensions Length × I.D. (mm)
ULTRON ES-OVM.G	5	10 × 4.0
ULTRON ES-OVM Prep.G (preparative guard column)	10	15 × 8.0
ULTRON ES-OVM (guard cartridges (x2))	5	5 × 2.0
Guard cartridge holder (with adaptor)	For 5 × 2.0 mm columns	
ULTRON ES-OVM (guard cartridges (x 2))	5	10 × 4.6
Guard cartridge holder (with adaptor)	For 10 × 4.6 mm columns	

### ● Protecting columns using guard columns

By using guard columns together with analytical columns, analytical columns can be protected from contamination, enabling longer column lifetimes. The use of guard columns is particularly recommended when analyzing samples containing large amounts of impurities.

Guard cartridges are also available when guard columns are changed frequently.

Please use guard columns with the same or smaller inner diameter than that of the analytical column being used.

### ● ULTRON ES-OVM-C (new product)

Clopidogrel bisulfate evaluation conditions based on USP standards employed for shipping inspections.

(Performance reports included with columns.)



Please feel free to contact us with questions related to analyses. Column screening services and demonstration columns are also available.

Please be aware that specifications and prices are subject to change without prior notification.