


SHINWA EXPRESS

Application

VOL. 215

Analysis of N₂O and NO₂ using Sunpak-A

Sunpak-A is an ethylvinylbenzene-divinylbenzene polymer packing material. Shown below is an analysis of N₂O and NO₂ performed using a column packed with a Sunpak-A packing material.

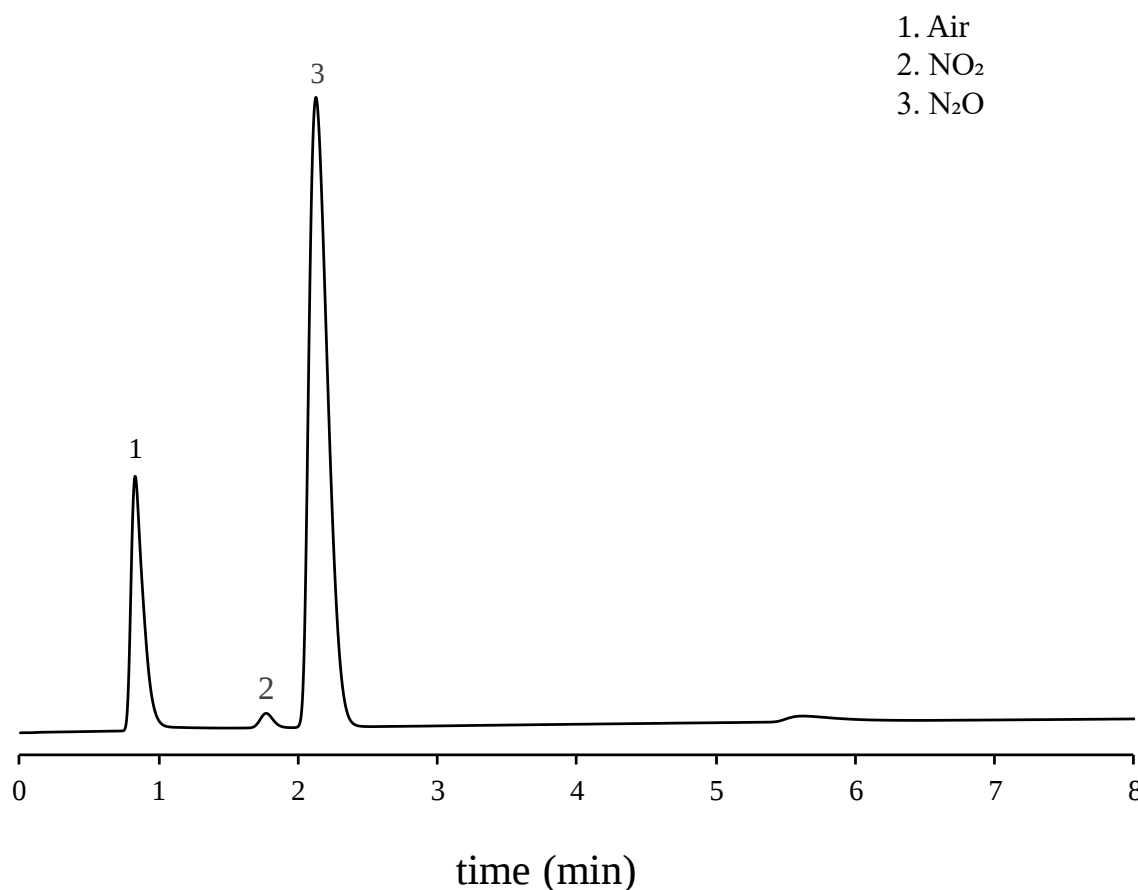

Column: Sunpak-A 80/100 mesh 2.0 m x 3.0 mm I.D. Glass
Column temp.: 50°C
Injection port temp.: 80°C

Detector: TCD 80°C
Carrier gas: 30 mL/min He
Sample vol.: 0.1 mL

**SHINWA CHEMICAL INDUSTRIES LTD.**

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>