

SHINWA EXPRESS


Application

VOL. 219

Analysis of Denatured alcohols using Sunpak-H

Sunpak-H is a styrene-divinylbenzene copolymer packing material.

Shown below is an analysis of denatured alcohols performed using a column packed with a Sunpak-H packing material.

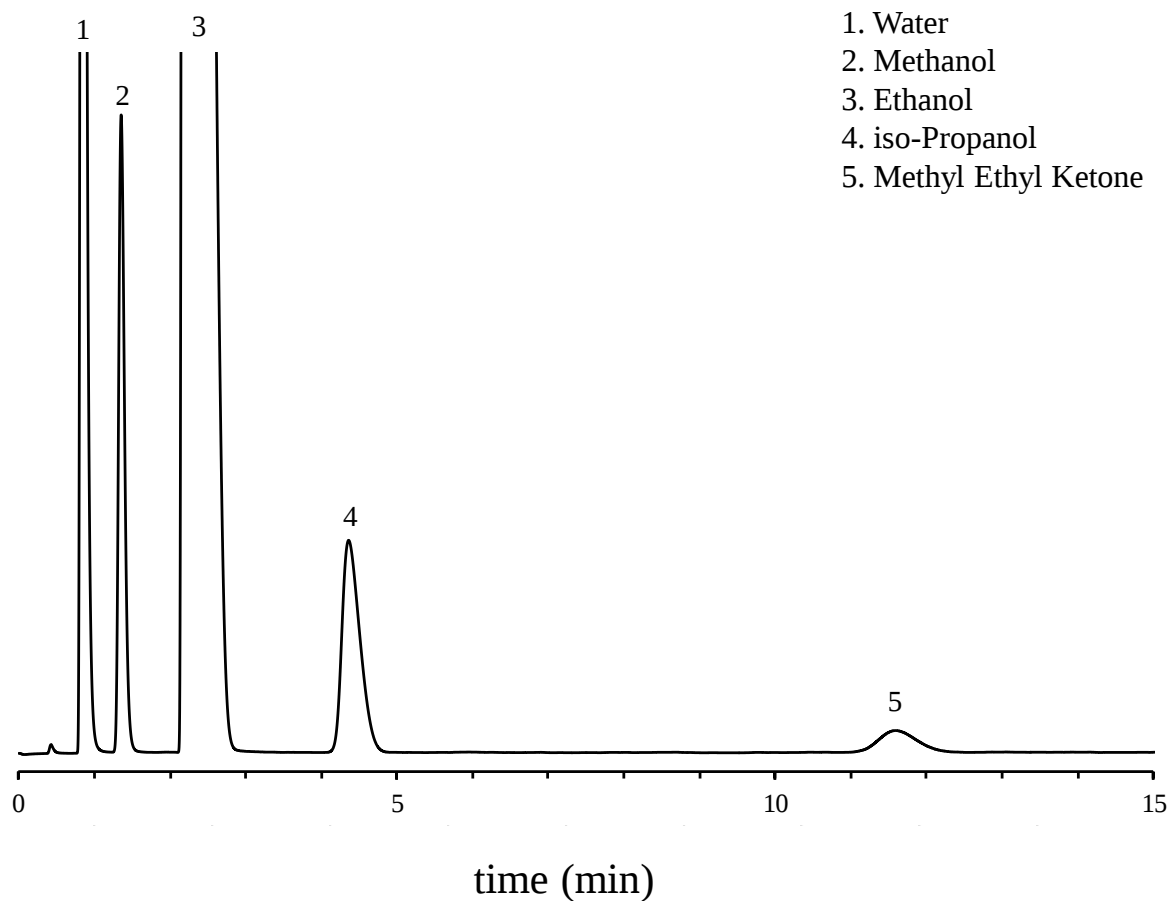

Column: Sunpak-H 60/80 mesh 2.0 m x 3.0 mm I.D. Glass

Column temp.: 120°C

Injection port temp.: 200°C

Detector: TCD 200°C

Carrier gas: 50 mL/min He

Sample vol.: 1.0 µL

**SHINWA CHEMICAL INDUSTRIES LTD.**

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan

TEL: +81-75-621-2360 FAX: +81-75-602-2660

E-mail: info@shinwa-cpc.co.jp

Website: <https://shinwa-cpc.co.jp/en/>