


SHINWA EXPRESS

Application

VOL. 223

Analysis of Inorganic gases using Shincarbon-ST

Shown below is an analysis of inorganic gases including propane and propylene performed using a column packed with a Shincarbon-ST packing material.

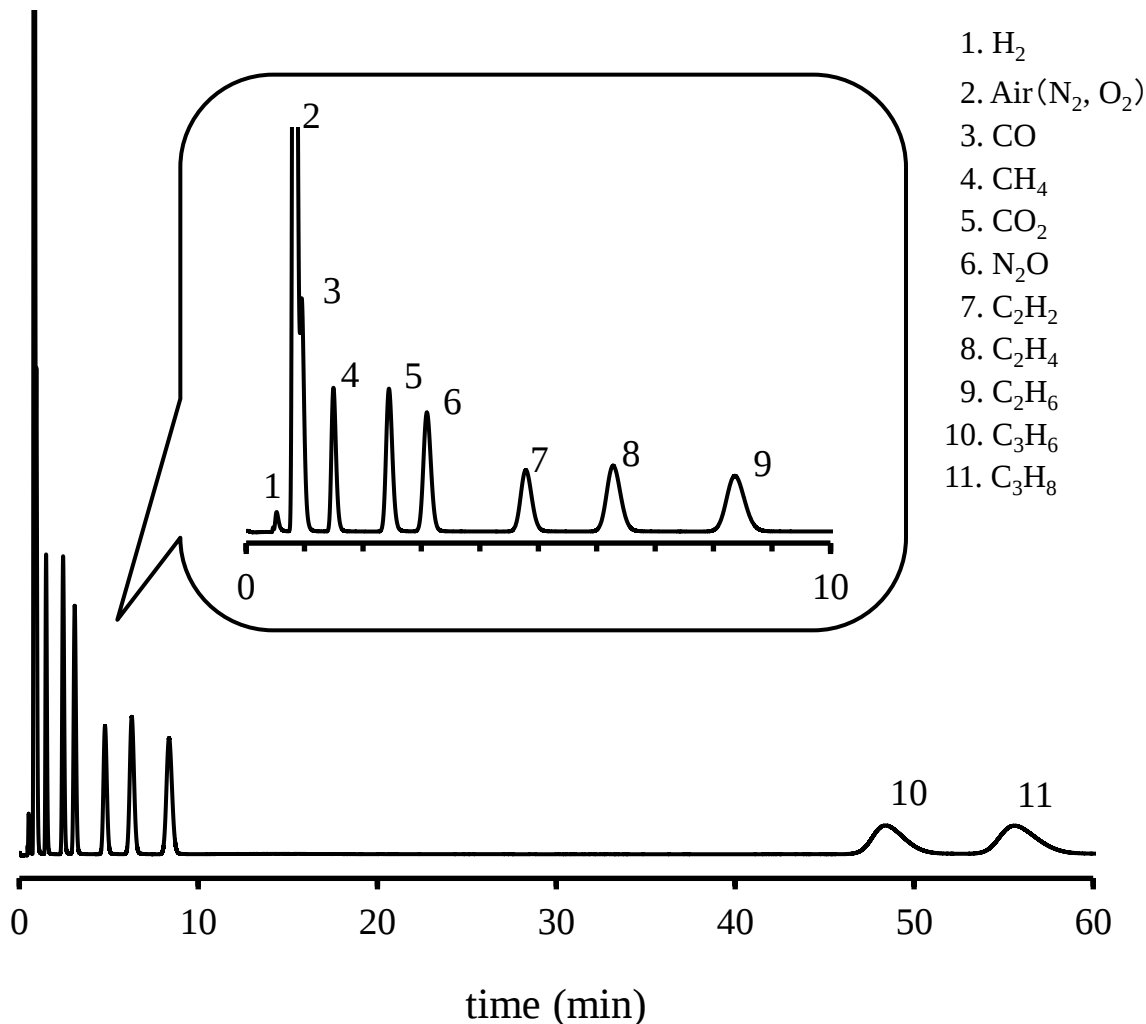

Column: Shincarbon-ST 50/80 mesh 2.0 m x 3.0 mm I.D. SUS

Detector: TCD 200°C

Column temp.: 150°C

Carrier gas: 50 mL/min He

Injection port temp.: 200°C

Sample vol.: 0.5 mL

**SHINWA CHEMICAL INDUSTRIES LTD.**

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>