

SHINWA EXPRESS


Application

VOL. 240

Analysis of Organic solvents (20 types) using Sunpak-E

Sunpak-E is a divinylbenzene-ethylvinylbenzene-ethylene glycol dimethacrylate copolymer polymer packing material.

Shown below is an analysis of organic solvents (20 types) performed using a column packed with a Sunpak-E packing material.

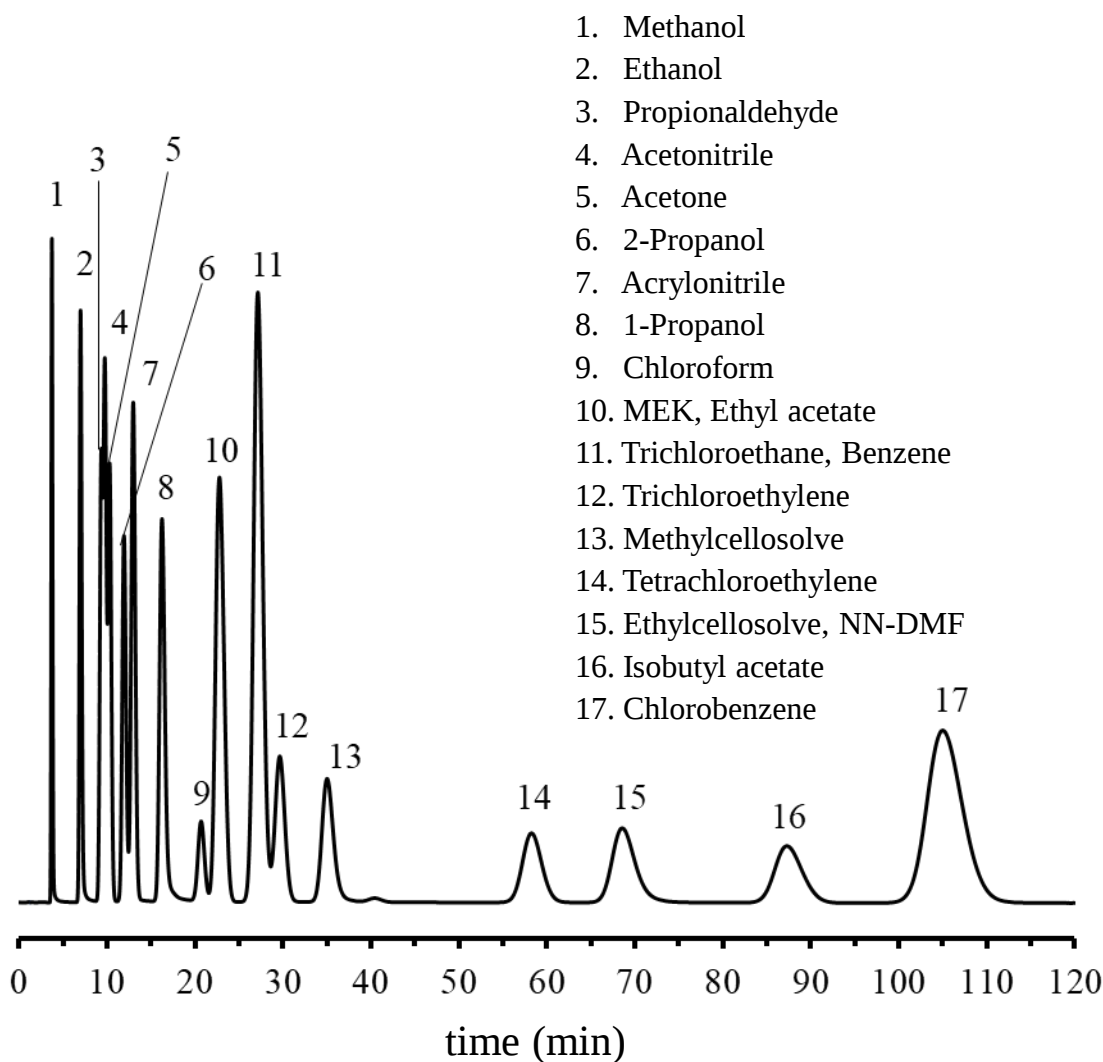

Column: Sunpak-E 80/100 mesh 2.0 m x 3.0 mm I.D. Glass

Column temp.: 175°C

Injection port temp.: 190°C

Detection: FID 190°C

Carrier gas: 20 mL/min N₂

Sample vol.: 0.2 µL



SHINWA CHEMICAL INDUSTRIES LTD.

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>