


SHINWA EXPRESS

Application

VOL. 244

Analysis of Organic solvents using Sunpak-A

Sunpak-A is an ethylvinylbenzene-divinylbenzene polymer packing material. Shown below is an analysis of organic solvents in water performed using a column packed with a Sunpak-A packing material.

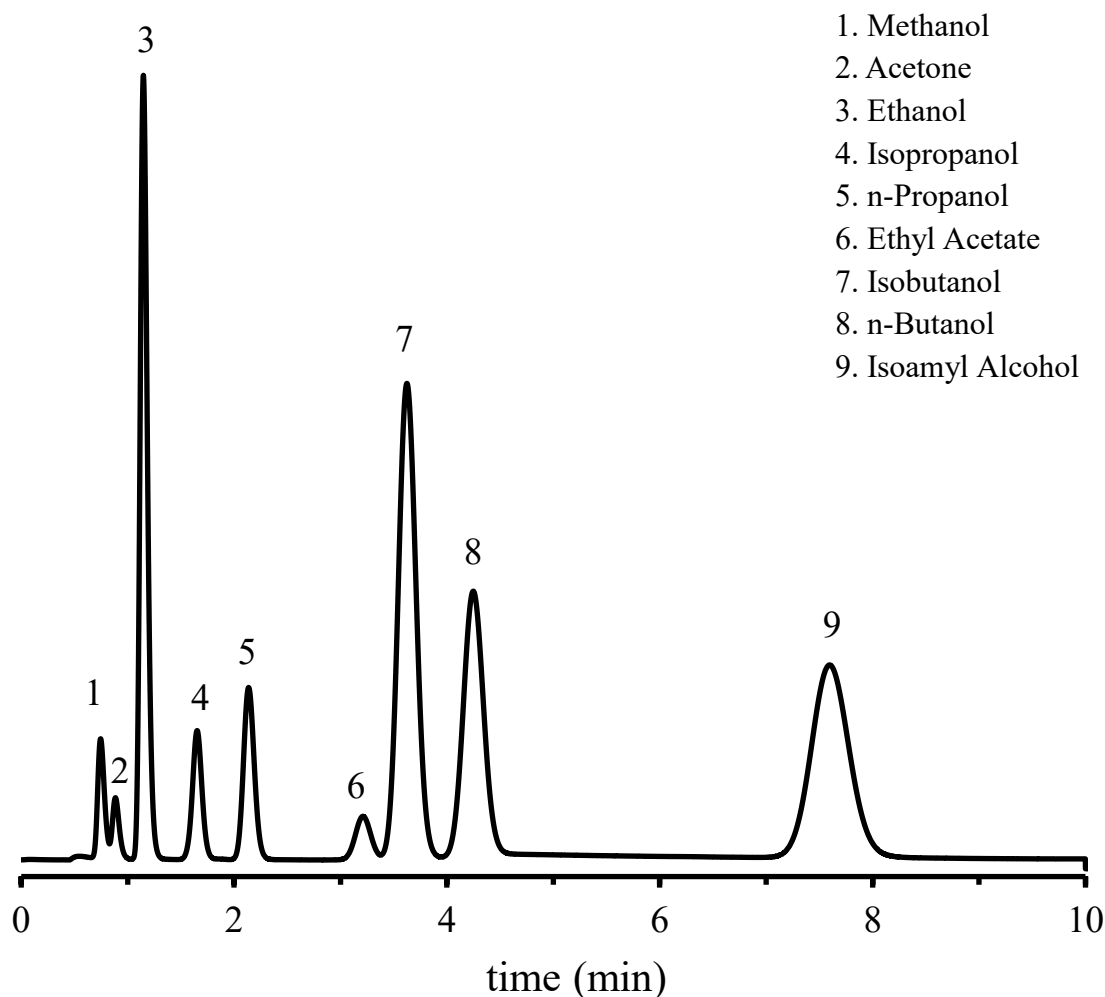

Column: 5% Thermon-1000 Sunpak-A 50/80 mesh
2.0 m x 3.0 mm I.D. Glass
Column temp.: 180°C
Injection port temp.: 250°C

Detection: FID 250°C
Carrier gas: 50 mL/min N₂
Sample vol.: 1.0 µL

**SHINWA CHEMICAL INDUSTRIES LTD.**

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>