

SHINWA EXPRESS


Application

VOL. 249

Analysis of Organic solvents (17 types) using Shinwasorb

Shinwasorb is a porous silica gas chromatography support.

Shown below is an analysis of organic solvents (17 types) performed using a column packed with a Shinwasorb packing material.

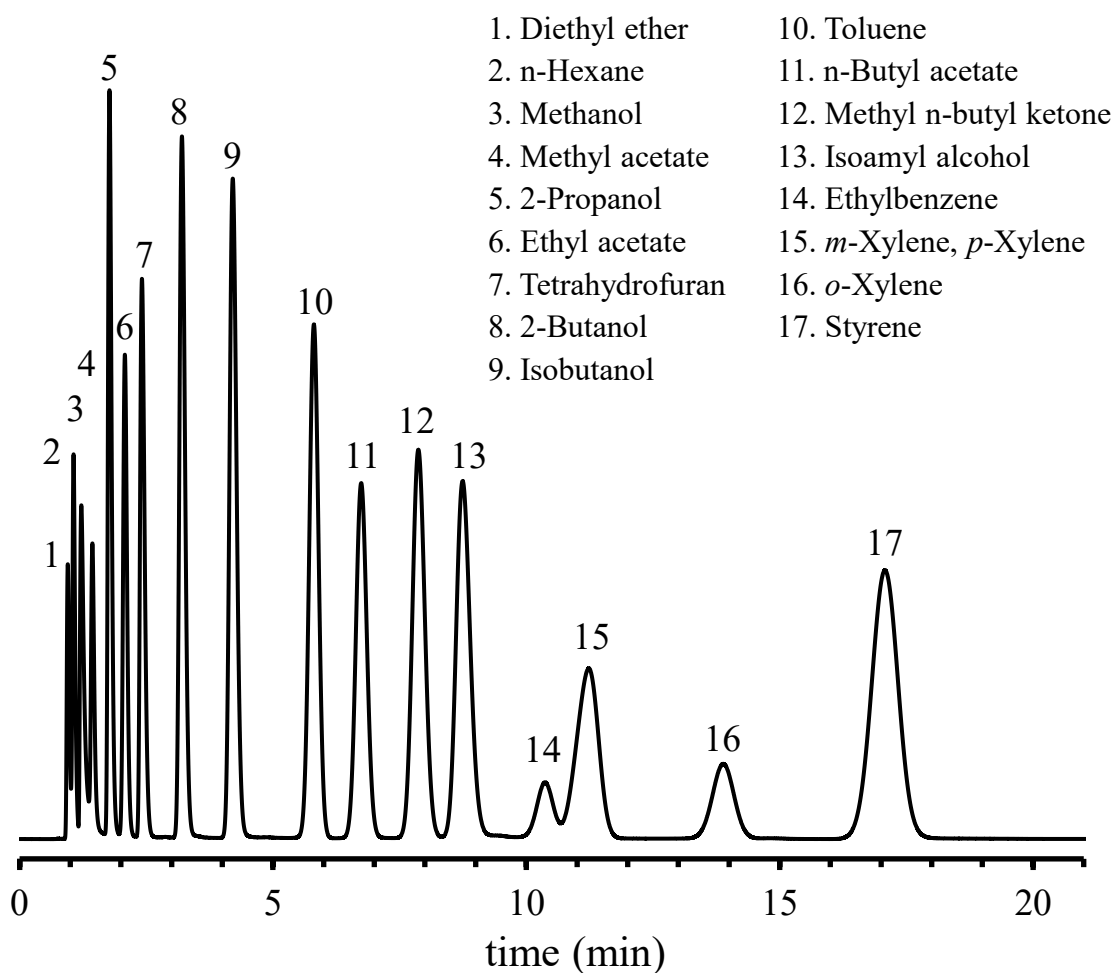

Column: 12% SBS-120 Shinwasorb-S 80/100 mesh
2.0 m x 3.0 mm I.D. Glass

Column temp.: 100°C

Injection port temp.: 200°C

Detection: FID 200°C

Carrier gas: 50 mL/min N₂

Sample vol.: 1.0 µL



SHINWA CHEMICAL INDUSTRIES LTD.

50-2 Kagekatsu-cho, Fushimi-ku
Kyoto 612-8307 Japan
TEL: +81-75-621-2360 FAX: +81-75-602-2660
E-mail: info@shinwa-cpc.co.jp
Website: <https://shinwa-cpc.co.jp/en/>